

FREQUENTLY ASKED QUESTIONS

HOW DOES IPL-SQ TREAT SKIN PIGMENTATION?

IPL-Sq's true square pulse system is delivered in waves at the right intensity to act on the darkened areas of pigment, minimising the risk to harm the skin and surrounding tissue.

WHAT RESULTS CAN I EXPECT?

You may experience some initial darkening of the treated areas. This is temporary, and after about 5-7 days, they will begin to fade. They will become progressively lighter with each treatment. Usually 3-4 treatments are required to reduce most of the pigment.

WILL IPL-SQ SUIT MY SKIN TYPE?

IPL-Sq can be used successfully on most skin types—please consult your medical practitioner for confirmation. As well as addressing pigmentation, IPL-Sq is an excellent choice for all-round skin rejuvenation treatment.

HOW WILL THE TREATMENT FEEL?

You will be aware of the pulses of light flashing. You will feel a slight sting, as if a rubber band is being snapped against your skin, but this is minimal, thanks to IPL-Sq's built-in adjustable contact cooling.

For best ongoing results, you should always protect your skin from direct sun exposure. Always cover up and use a good SPF sunblock when out in the sun.

IPL-SQ

FADE SUN DAMAGE & SKIN PIGMENTATION

TRUE SQUARE PULSE SYSTEM

- Targets darkened areas of pigmentation
- Breaks up areas of pigmentation, allowing the body to remove it naturally
- Surrounding tissue is not damaged
- Most skin types can be treated
- No medication necessary

Clinic details:

UA ULTRA
AESTHETICS

ultraa.com.au

Ultra Aesthetics is a division
of the Cryomed Group.

Cryomed
Aesthetics

cryomed.com.au | 1300 346 448

Please consult your medical practitioner to see if you are suitable for treatment. Potential risks and side effects need to be discussed with your medical practitioner before undergoing treatment with the IPL-Sq.

© Ultra Aesthetics 2016

CLEAR AWAY PIGMENTATION



**RESTORE YOUR
CLEAR, HEALTHY
COMPLEXION**

ARE YOU TIRED OF BLOTCHY SKIN?

IPL-Sq targets and destroys the pigmentation in freckles, pigmented lesions, sunspots and other dark pigmentation on your skin.

STEP 1

The skin is cleaned and a clear coupling gel is applied to the treatment area. Your eyes will be covered and you will be instructed to keep them closed.

STEP 2

Light energy is delivered to the skin with the most advanced square-pulse technology. The operator will move the head over every area of pigmentation in the treatment area. You will sense the flash of light and feel a sensation like a rubber band being snapped against your skin.

STEP 3

This causes the pigmented areas to heat up, then break down and flake off. The process is called photothermolysis.

In some circumstances, IPL-Sq treatment might need to be postponed. You should seek medical advice if any of the following apply to you:

- You have tanned skin
- You have a disease that is stimulated by light (e.g., Epilepsy)
- You are using photosensitive drugs (e.g., Roacutane)
- You are pregnant

THE IPL-SQ RESTORES YOUR CLEAR COMPLEXION

A beautiful creamy complexion gives you a more youthful appearance. Your skin gets its colour from a pigment called melanin, which is produced by specialised skin cells.

If these cells become damaged or unhealthy, they can become highly erratic in their melanin production, leading to darker or lighter patches of skin. The dark patches are sometimes referred to as 'age spots', but younger people can also develop them.

Most common types of benign pigmented lesions can be treated with the IPL-Sq, with excellent efficacy and clinically proven results. IPL-Sq is advanced and acclaimed technology for non-invasive fading of skin pigmentation on the face and body.

IPL-Sq's bursts of intense pulsed light target the darkened areas of pigment in sun damage, causing it to heat up and then break down. This process is known as photothermolysis. As the dark areas break down, they flake off.

IPL-Sq treatment is fast and comfortable—with no downtime. You can return to your normal activities immediately after treatment, with predictable and dramatic results.

RESULTS YOU CAN SEE AND FEEL



Before




After



Before



After

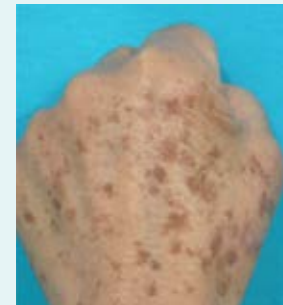

Before



After

* Individual results may vary based on the nature and severity of the skin condition and aggressiveness of treatment.