Ulfit

New transformed Face lifting



Body contouring System



What is Ulfit?


A non-invasive Micro Focused Ultrasound & Macro Focused Circular Ultrasound

FACE lifting and Body Contouring System





Trend & New Paradigm

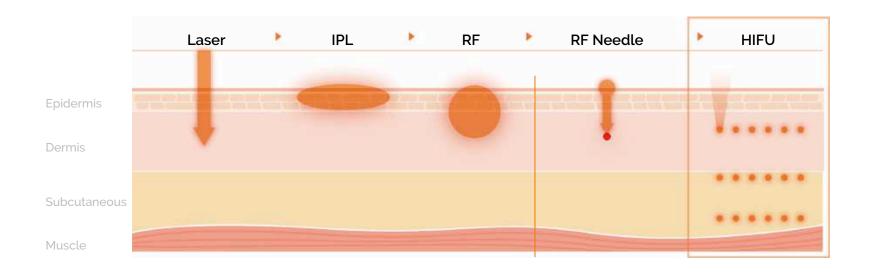

- Long downtime.

- Skin Damage / Hyper pigmentation 1st Co2 laser by Kumar Patel& Bell Labs 1964


- Short downtime.
- Less pain.
- Bleeding
- Efficacy on epide rmal layer
- No downtime.
- No skin surface dama ge at all
- Intensity focused ene rgy emitting on the ta rget area only.
- HIFU creates 65~75 d egrees of thermal co agulations in the skin depth of 1.5mm, 3mm , 4.5mm
- Without damaging sk in surface.



Mechanism

Micro Focused Ultrasound


Micro focused face cartridges deliver its energy precisely to the targeted area much faster than any others.

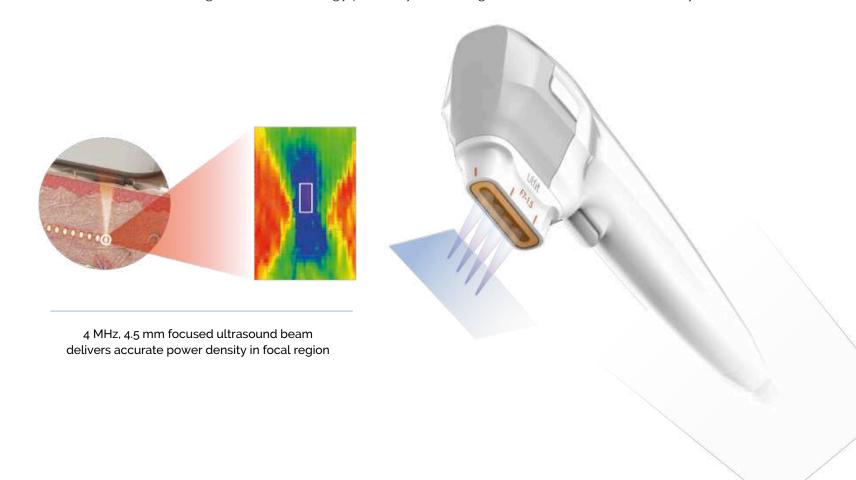

Mechanism

Macro Focused Circular Ultrasound

Macro Focused Circular body cartridge continuously delivers its ultrasound energy in circular way and it maximizes its treatment result in a short time. Motional movement on the treatment area also dramatically reduces pain level.

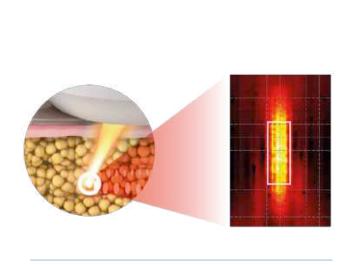

2 MHz, 13.0 mm focused ultrasound beam delivers higher power density in focal region





Face Cartridges

MFU (Micro Focused Ultrasound)

F7-1.5 mm

Non-Surgical Face Lifting

Face TX

F7-3.0 mm

F4-4.5 mm



Body Cartridges

MFCU (Macro Focused Circular Ultrasound)

Motional Body Contouring

Body TX






B2-13.0 mm



Face & Body Cartridges

MFU (Micro Focused Ultrasound) & MFCU (Macro Focused Circular Ultrasound)

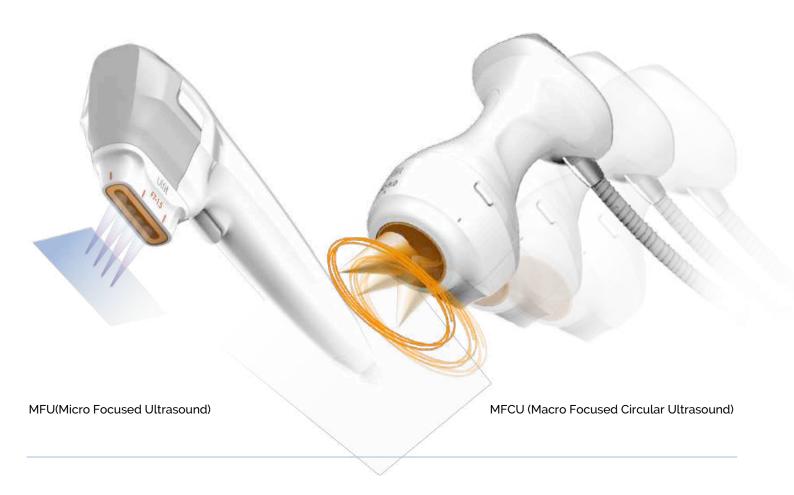




01

Fast shot speed of Face applicator

_

Higher Peak Power by strong Dual Engine

- Less pain
- Long and stable shot quality
- Better efficacy
- Higher max power(1.5 J/Face, 3.0J/Body)

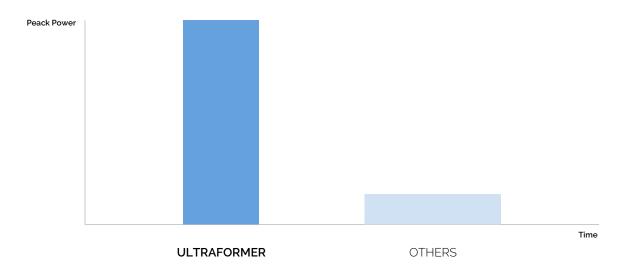



0.2

The world First Large Circular Motion Applicator

_

Creative concept of motional body applicator sprays out its ultrasound energy into targeted area by circular movement.






03

Extensive and Fast Energy Delivery

_

Specially designed transducer delivers numerous numbers of energy fastly in order to reduce treatment time.

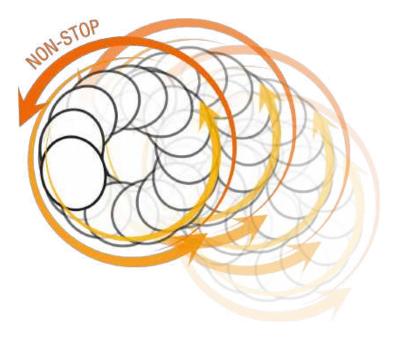



04

Extensive and Fast Energy Delivery

_

Higher Peak Power by strong Dual Engine

- Less pain
- Long and stable shot quality
- Better efficacy
- Higher max power(1.5 J/Face, 3.0J/Body)





05

Automatic Movable User Interface

_

Newly designed user interface automatically moves to right or left positions when operator holds the handpieces.

Automatic sensor system also recognizes the handpieces which is connected in each side.



Indication for use

Face TX

.



Body TX

_





Specification



Energy	MFU & MFCU
Fluence	0.1~3.0J (Face: up to 1.5J)
Cartridge	7Mhz, 4Mhz 2Mhz (1.5 mm, 3 mm, 4.5 mm) (13 mm)
Spacing	1.0 - 2.0 mm (0.5mm/step)
Length	5-25 mm (5 mm/step)
Electrical Requirements	AC100~240V, 50/60Hz
Weight	14.4 Kg
Dimensions	360 x 460 x 150 mm (WxDxH)



Treatment Procedure

Ulfit treatment is

Smarter & Faster





Total treatment for face and neck takes 15 – 20 minutes & body takes 20 – 30 minutes



Ulfit


Clinical Cases

Face Treatment






Body Treatment

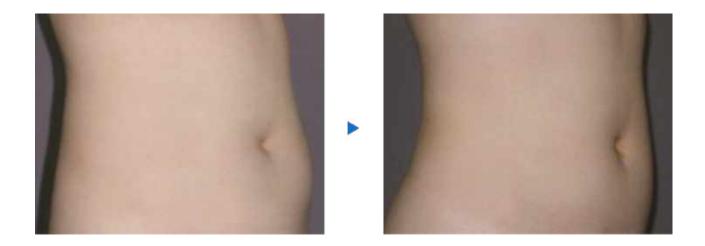



Body Treatment

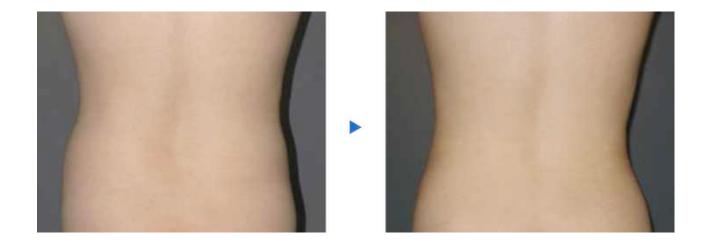



Body Treatment





Thank you!