

# Quick Operation Guide ULTRAFORMER MPT

ver 1.2



## 1. Press down the power button

Ensure the Main Power Switch is in the ON position. Press down the key button to turn the System ON.



## 2. Set the parameters

Set the parameters (Energy, Pitch, Length, Mode) by adjusting the buttons on the Graphical User Interface (GUI).



## 3. Apply ultrasound gel

Apply a thin layer of ultrasound transmission gel to the treatment area and cartridge. Use ultrasound transmission gel only, as other lubricants or lotions may damage the Cartridges or Hand-piece.



## 4. Apply the Cartridge

Place the Cartridge window on the treatment areas.




## 5. Commence treatment

Select the "STANDBY" button on the GUI. Once selected it will change to "READY". And pull the trigger (operating source: hands or foot switch) to emit ultrasound energy.



## 6. Check skin reactions and sensitivity of patient

After a few shots, check for skin reactions and sensitivity level of the patient. Proceed the treatment if there are no adverse events.

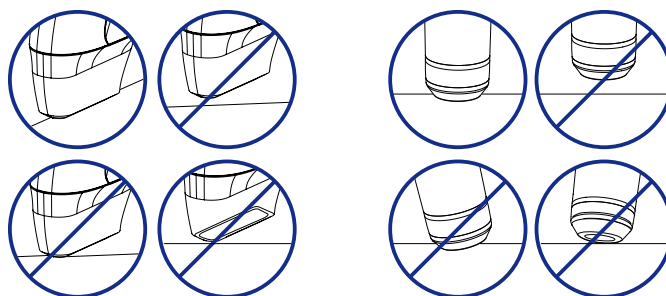

DO NOT apply shots without making proper contact between the cartridge and skin.

DO NOT press the shot button of the hand-piece without applying ultrasound gel onto treatment area.

## ○ Contraindication

The ULTRAFORMER MPT is contraindicated for use in patients with:

- Pregnancy or when breastfeeding
- Open wounds or lesions on the treatment area
- Severe or cystic acne on the treatment area
- Metal stents/electrical implants in the treatment area
- Bio-absorbable mechanical implants
- Skin infections of any type
- Hemorrhagic disorders or dysfunctions

## ○ Precautions for Patient

Determine whether treatment is appropriate for patients with the following conditions:

- Anticoagulant therapy
- Skin conditions caused by autoimmune diseases, skin cancer, and Herpes Simplex
- Diabetes or epilepsy
- Severe skin diseases
- Decortication or heart disease
- Hypersensitiveness or hypertrophic scars
- Hemostatic dysfunctions

## ○ Precautions for Procedure

- Do not use the ULTRAFORMER MPT without completion of training by an authorized Classys representative
- To prevent unauthorized use of the device, ensure access is limited strictly to trained and certified personnel only
- To avoid possible injury during treatment, do not apply the ULTRAFORMER MPT hand-pieces near the patient's eyes
- Avoid applying treatment on the thyroid gland, thyroid cartilage, and trachea
- Avoid applying treatment directly on the breast tissue and major vessels
- If any malfunctions or adverse events occur, immediately stop the treatment procedure
- Always apply ultrasound gel on target areas before proceeding with the ULTRAFORMER MPT treatment
- Before conducting shots, check that the cartridge attached on the transducer is applied on the skin properly
- Immediately stop treatment if the surface of the skin burns or the outer surface layer of skin rises by swelling

## ○ Precautions for System Safety

- To avoid risk of electric shock, always inspect the cartridge, hand-piece and cable before use.
- Do not use a cable or cartridge that has been damaged or is leaking fluid.
- Do not remove the covers attached to the hand-piece.
- The ULTRAFORMER MPT contains no user serviceable components. If the system requires service, contact a local distributor or Classys representative. Modification of this equipment is not allowed.
- Do not lean over the equipment obliquely and do not give physical external shock to the equipment as it may cause permanent system damage.
- Do not use the equipment with anesthetic gases or flammable materials in the surrounding areas.




For additional information refer to the Operation Manual of the ULTRAFORMER MPT


# Quick Operation Guide ULTRAFORMER MPT

ver 1.2

## Graphical User Interface

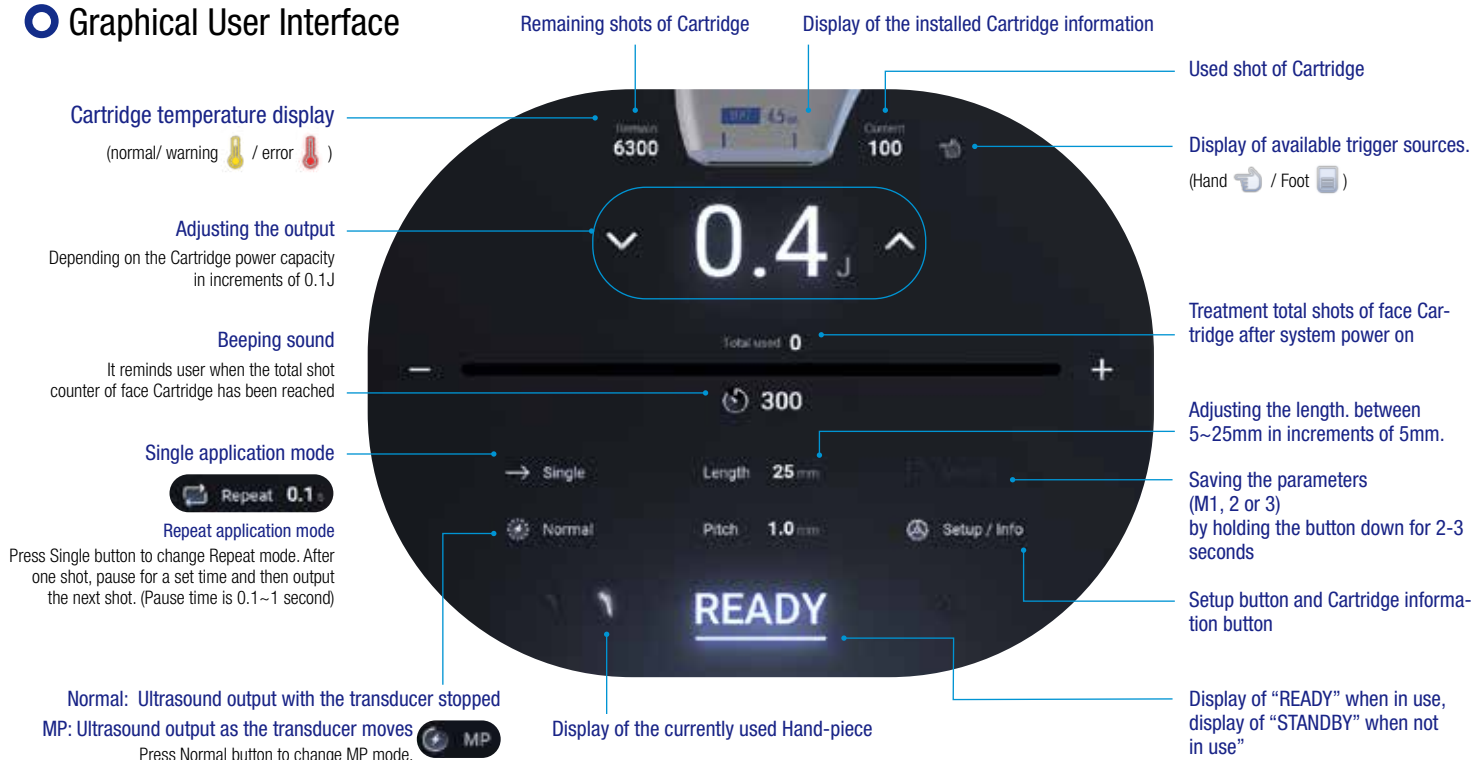

## Fine Lines & Wrinkles

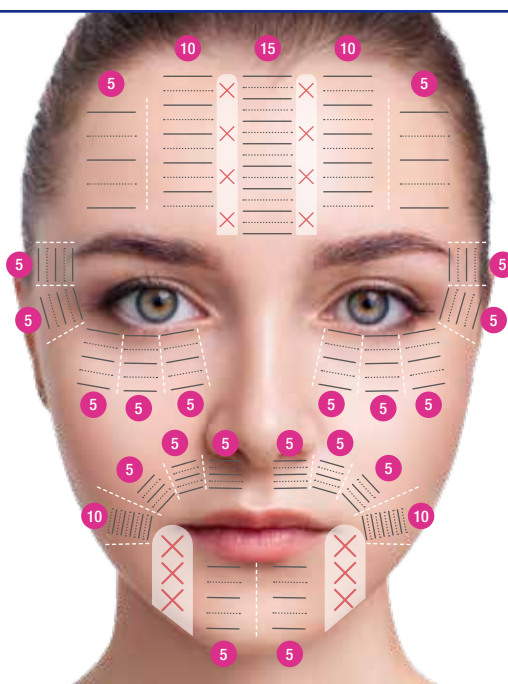


Area	Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
Forehead	Normal	2.0mm	0.1-0.3J	1.0-1.5mm	1	45
	MP		0.1-0.3J	-	1	45
Around the Eyes	Normal		0.1-0.2J	1.0-1.5mm	1	50
	MP		0.1-0.2J	-	1	50
Around the Mouth	Normal		0.2-0.3J	1.0-1.5mm	1	60
	MP		0.2-0.3J	-	1	60

- Do not apply treatment directly on the supraorbital or temporal nerves.
- Adjust the shot length between 10mm and 15mm as it is difficult for the transducer to make firm contact with both the periorbital region and forehead.
- Have patients close their eyes when treating the periorbital region.
- When treating near the eyebrow line, pull the transducer tighter towards the forehead to reduce pain.
- For treatment applications with 1 pass, apply the treatment pattern as directed.

# Quick Operation Guide ULTRAFORMER MPT

ver 1.2



## Forehead Wrinkles

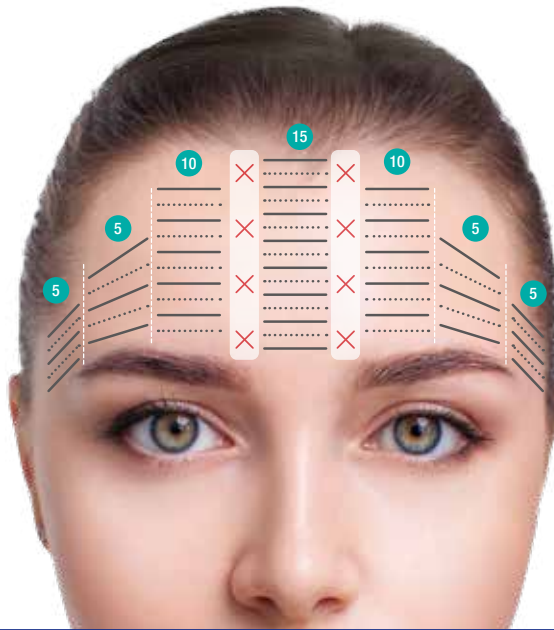

Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
Normal	1.5mm	0.1-0.2J	1.0-1.5mm	1	55
MP		0.1-0.2J	-	1	55

- Do not apply treatment directly on the supraorbital or temporal nerves.
- Adjust the shot length between 10mm and 15mm as it is difficult for the transducer to make firm contact with both the periorbital region and forehead.
- Have patients close their eyes when treating the periorbital region.
- When treating near the brow line, pull the transducer tighter towards the forehead to reduce pain.
- For treatment applications with 1 pass, apply the treatment pattern as directed.

## Lower Face Contouring / Marionette Line

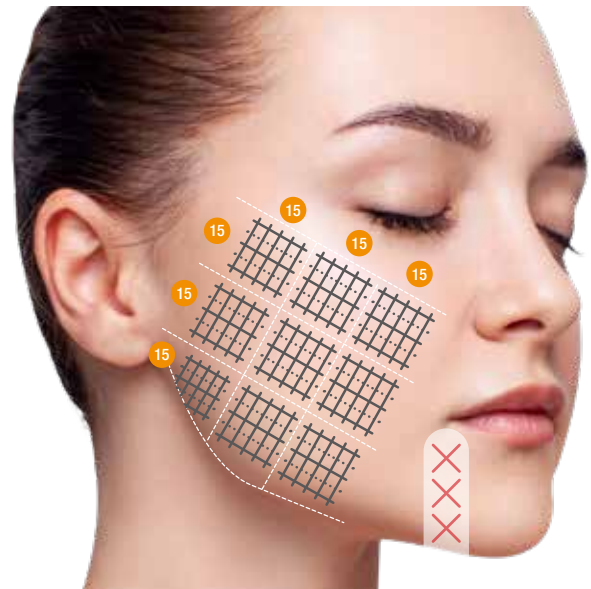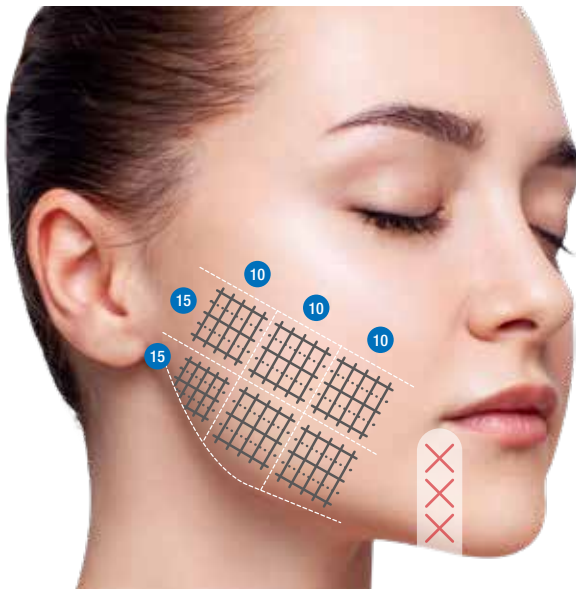



Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
Normal	4.5mm	0.5-0.6J	1.0-1.5mm	1	60
MP		0.4-0.5J	-	1	60
Normal	3.0mm	0.4-0.5J	1.0-1.5mm	1	90
MP		0.4-0.5J	-	1	90

- Do not apply treatment directly on the marginal mandibular nerve.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 4.5mm ► 3.0mm
- For treatment applications with 2 passes, practitioners are recommended to apply treatment in both horizontal and vertical patterns.
- For treatment applications with 1 pass, apply the treatment pattern as directed.

# Quick Operation Guide ULTRAFORMER MPT

ver 1.2

## ● Double Chin & Neck

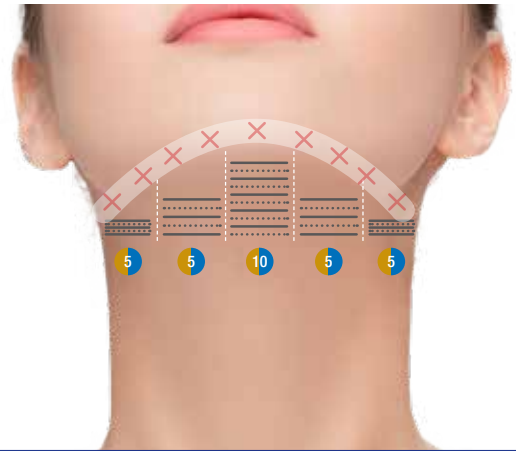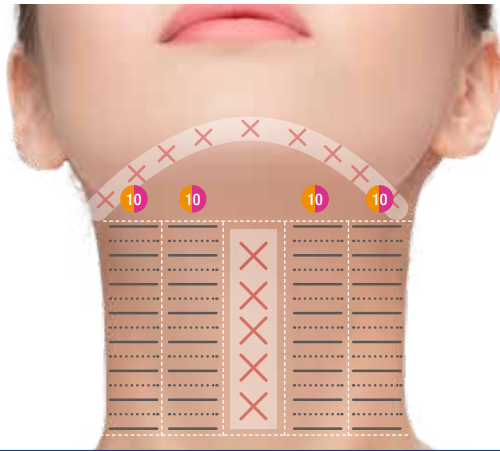

Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
Normal	3.0mm	0.4-0.5J	1.0-1.5mm	2	50
MP		0.3-0.4J	-	2	50
Normal	2.0mm	0.1-0.3J	1.0-1.5mm	2	100
MP		0.1-0.3J	-	2	100

- Do not apply treatment directly on the marginal mandibular nerve and thyroid cartilage.
- Do not apply excess pressure on the neck to prevent discomfort during procedures.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 3.0mm ▶ 2.0mm
- For treatment applications with 1 pass, apply the treatment pattern as directed.

Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
Normal	6.0mm	0.6-0.8J	1.0-1.5mm	2	50
MP		0.6-0.8J	-	2	50
Normal	4.5mm	0.4-0.6J	1.0-1.5mm	2	60
MP		0.4-0.6J	-	2	60

- Do not apply treatment directly on the marginal mandibular nerve and thyroid cartilage.
- Do not apply excess pressure on the neck to prevent discomfort during procedures.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 6.0mm ▶ 4.5mm
- For treatment applications with 1 pass, apply the treatment pattern as directed.

## ● Glowing and SMAS Lifting (using Booster hand-piece)



Mode	Cartridge(BOOSTER)	Energy	Pitch	Shots
Normal	1.5mm	0.2-0.3J	1.0-1.5mm	150-200
MP		0.2-0.3J	-	150-200
Normal	3.0mm	0.3-0.5J	1.0-1.5mm	100-150
MP		0.3-0.5J	-	100-150
Normal	4.5mm	0.4-0.5J	1.0-1.5mm	100-150
MP		0.4-0.5J	-	100-150

- Do not apply treatment directly on the marginal mandibular nerve.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 4.5mm ▶ 3.0mm ▶ 1.5mm
- For treatment applications with 1 pass, apply the treatment pattern as directed.

# Quick Operation Guide ULTRAFORMER MPT

ver 1.2


## Upper Arms



Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
MP	6.0mm	0.6-0.8J	-	2	100
MP	4.5mm	0.6-1.0J	-	2	100
MP	3.0mm	0.4-0.6J	-	2	100

- Do not apply treatment directly on the ulnar nerve.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 6.0mm ► 4.5mm ► 3.0mm
- For treatment applications with 2 passes, practitioners are recommended to apply treatment in both horizontal and vertical patterns.
- This parameter is applied for each arm.

## Axilla

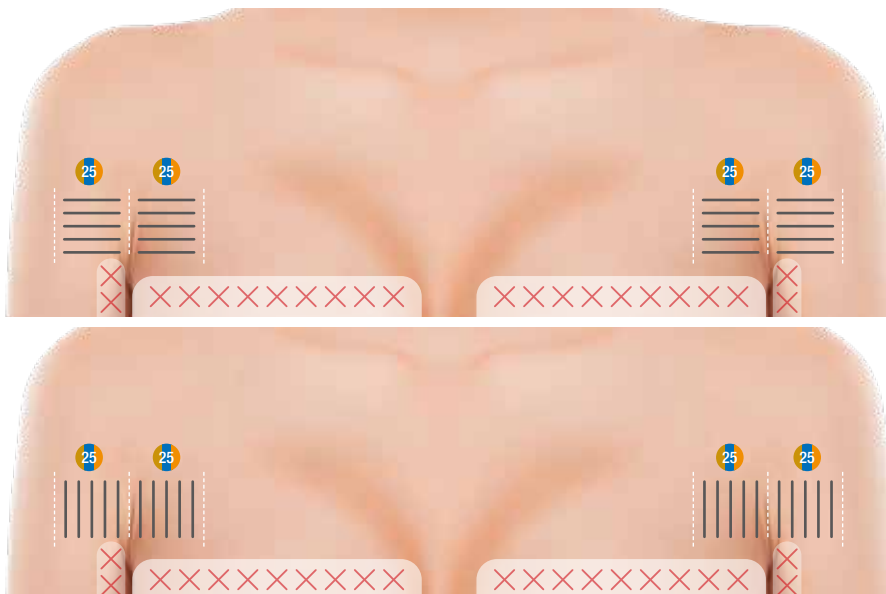

Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
MP	6.0mm	0.3-0.5J	-	2	100
MP	4.5mm	0.3-0.5J	-	2	100
MP	3.0mm	0.3-0.5J	-	2	100


- Do not apply treatment directly on breast tissue or the ulnar nerve.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 6.0mm ► 4.5mm ► 3.0mm
- For treatment applications with 2 passes, practitioners are recommended to apply treatment in both horizontal and vertical patterns.
- This parameter is applied for each axilla.



# Quick Operation Guide ULTRAFORMER MPT

ver 1.2


## Decollete

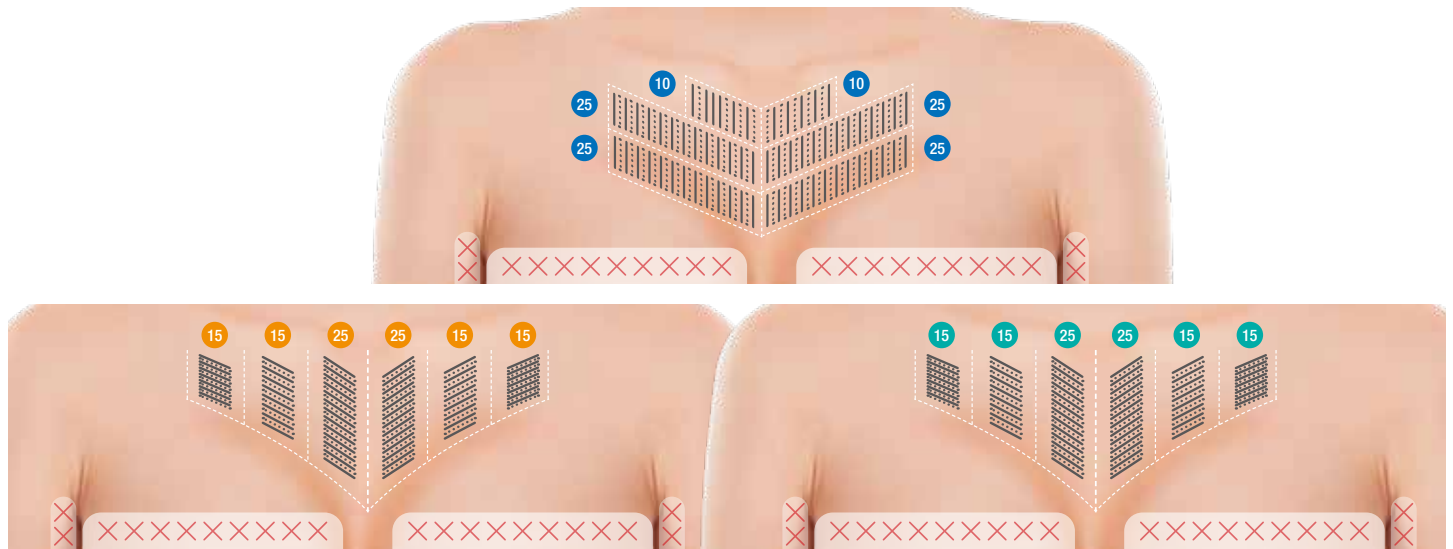

Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
Normal + MP	4.5mm	0.3-0.5J	1.0-1.5mm	2	120
Normal + MP	3.0mm	0.3-0.5J	1.0-1.5mm	1	100
Normal + MP	1.5mm	0.1-0.2J	1.0-1.5mm	1	100

- Do not apply treatment directly on the clavicle, breast tissue, or breast implants.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 4.5mm ▶ 3.0mm ▶ 1.5mm
- For treatment applications with 1 pass, apply the treatment pattern as directed.
- 4.5mm Cartridge is not recommended for the patient with thin skin layer.

## Abdomen

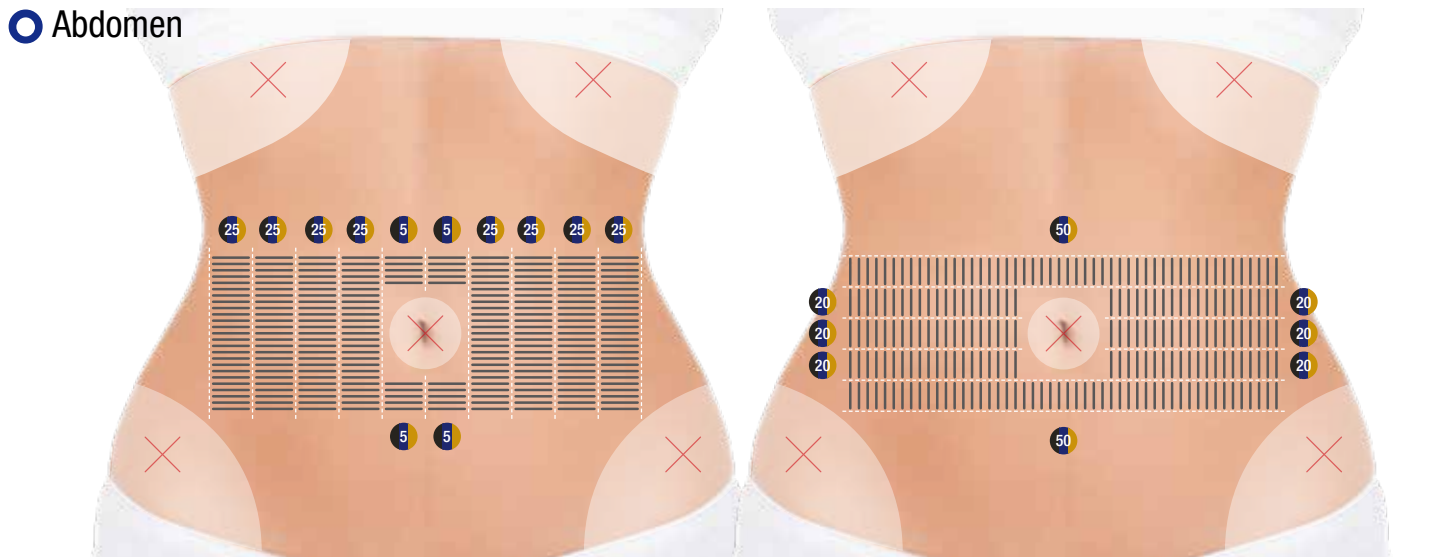


Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
MP	13.0mm	1.2-1.6J	-	2	300
MP	9.0mm	1.2-1.6J	-	2	300
MP	6.0mm	0.8J-1.2J	-	2	300

- Do not apply treatment directly on the navel, ribs, or pelvis.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 13.0mm ▶ 9.0mm ▶ 6.0mm
- For treatment applications with each pass, practitioners are recommended to apply treatment with both horizontal and vertical treatment patterns as shown above. Practitioners may choose to treat between either the horizontal or vertical pattern to be applied twice.

# Quick Operation Guide ULTRAFORMER MPT

ver 1.2


## Love Handles

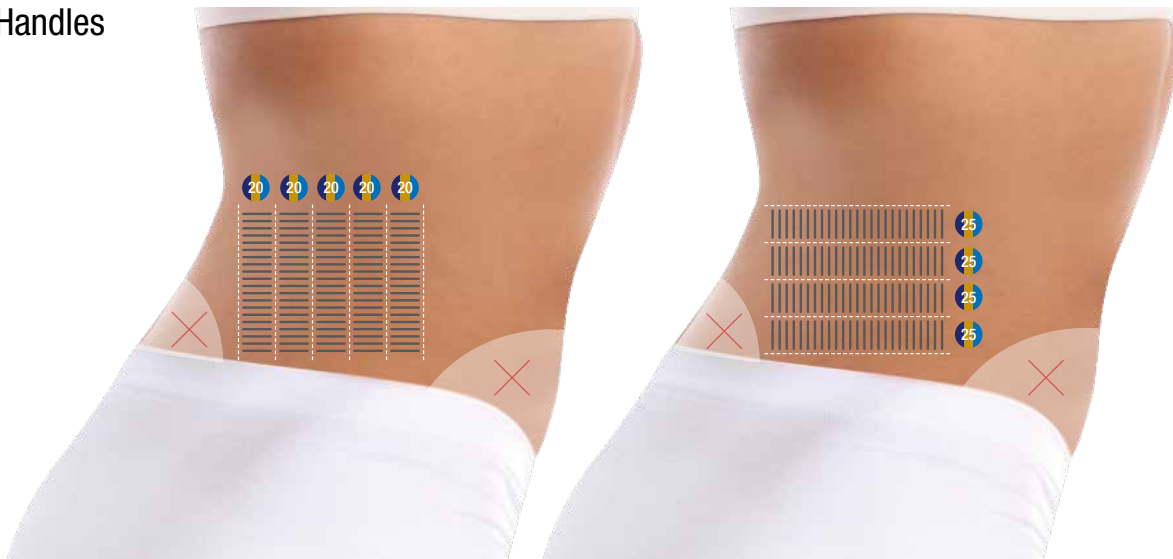

Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
MP ○	9.0mm	0.8-1.5J	-	2	200
MP ○	6.0mm	0.8-1.2J	-	2	200
MP ○	4.5mm	0.8-1.0J	-	2	200

- Do not apply treatment directly on the posterior superior iliac spine (PSIS) or pelvis.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 9.0mm ▶ 6.0mm ▶ 4.5mm
- For treatment applications with each pass, practitioners are recommended to apply treatment with both horizontal and vertical treatment patterns as shown above. Practitioners may choose to treat between either the horizontal or vertical pattern to be applied twice.

## Knees

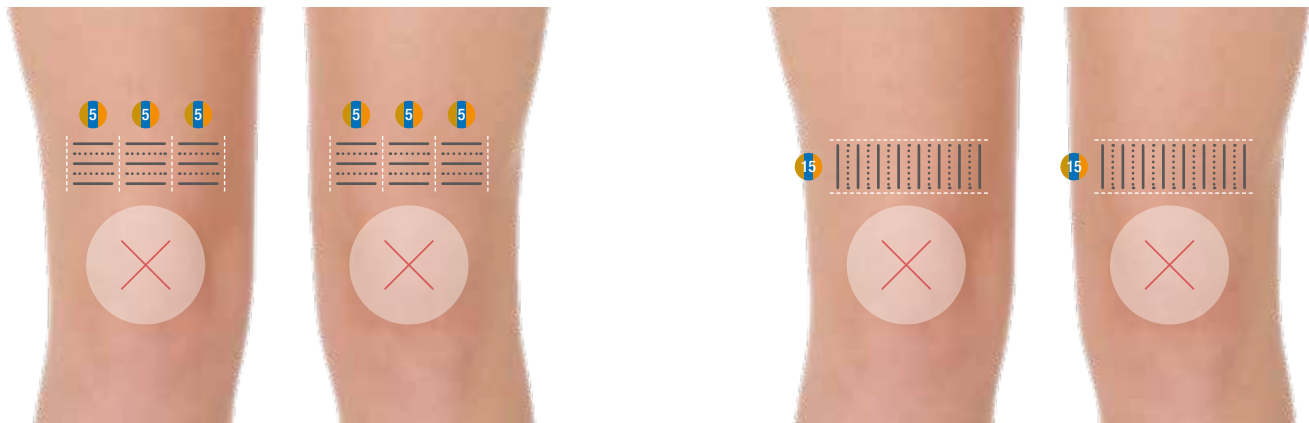


Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
Normal ○ + MP ○	6.0mm	0.6-0.8J	1.0-1.5mm	1	50
Normal ○ + MP ○	4.5mm	0.6-0.8J	1.0-1.5mm	1	50
Normal ○ + MP ○	3.0mm	0.5-0.6J	1.0-1.5mm	1	50





- Do not apply treatment directly on the kneecap or nerve regions.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 6.0mm ▶ 4.5mm ▶ 3.0mm
- For treatment applications with each pass, practitioners are recommended to apply treatment in both horizontal and vertical patterns.

# Quick Operation Guide ULTRAFORMER MPT

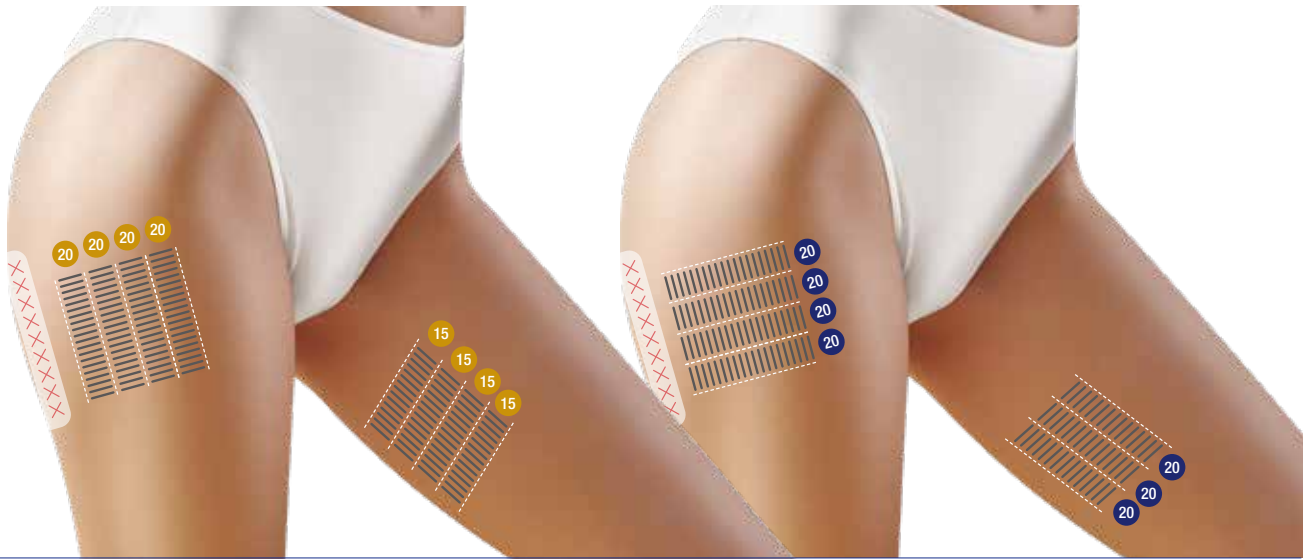ver 1.2

## Thighs



Mode	Cartridge(MPT)	Energy	Pitch	Passes	Shots
MP 	 9.0mm	0.6-0.8J	-	2	280
MP 	 6.0mm	0.6-0.8J	-	2	280

- Do not apply treatment directly on the sciatic nerve.
- When treating multiple layers, treat from deep to shallow in sequential order.  
Ex: 9.0mm ► 6.0mm
- For treatment applications with 3 passes, practitioners are recommended to apply treatment with both horizontal and vertical treatment patterns as shown above. Practitioners may choose to treat between either the horizontal or vertical pattern to be applied twice.