
Secret RF cases

(with treatment parameters)

SECRET: Acne scars

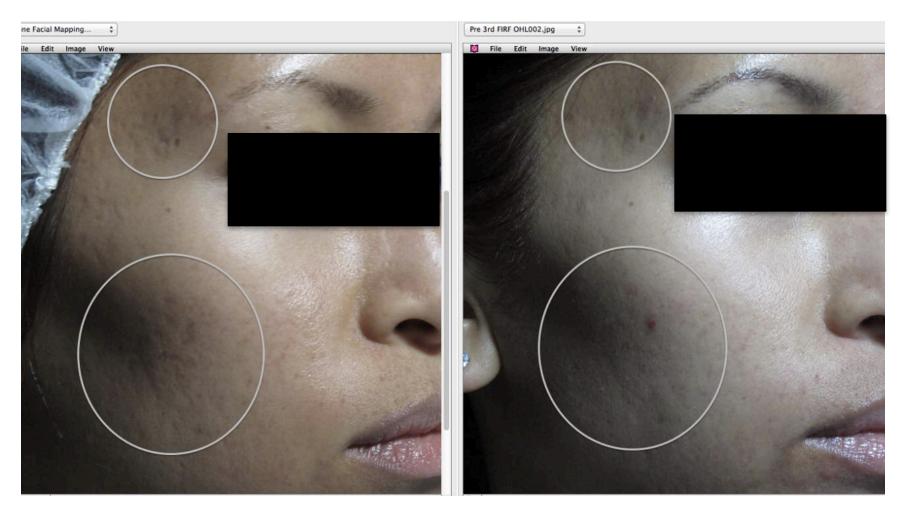

baseline 1mp2

SECRET: Acne scars


- Presenting Problem: Acne
- **1. History:** cheeks 60%/250/3.5/500 3x passes
- temples 60%/250/2.0/500 3x passes
- nose 60%/250/2.0/500 2x passes
- full face 60/250/3.0/500
- pt tolerated well with Siberian cooler
- stated before treatment she saw little improvment from last treatment, settings increased

- 2. History: 1st FIRF glabella, temples and cheeks
- 50%/250/3.0/500 Temles and cheeks
- 50%/200/3.0/400 Glabella

Well toletated with Siberian cooler


pt wants stretch marks to treat.

SECRET - Acne scars

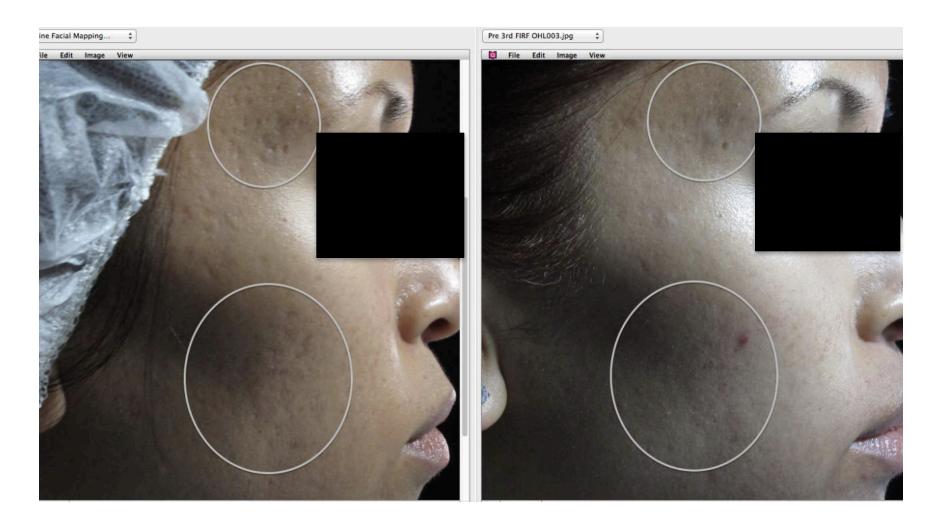

baseline 1mp2

SECRET - Acne scars

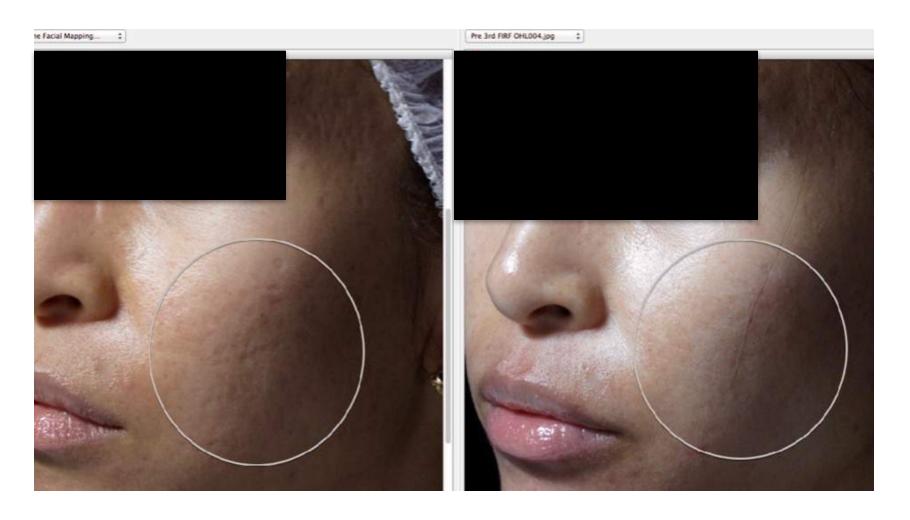

baseline 1mp2

SECRET - Acne scars



Laugh lines / crepiness



baseline 6wp2

Laugh lines / crepiness

- Presenting Problem: Laugh lines / crepiness
- 1. History: 2nd FIRF on Lower face
- then FIRF 60%/200/2.0/400 1-2 passes
- 3 passes on lines
- 2. History: 1st FIRF on Lower face
- then FIRF 50%/200/2.0/400 2 passes
- next treatment in a month
- Siberian cooler
- CA for 45min


Medial cheek hollowing



baseline 3mp3

- Presenting Problem: Medial cheek hollowing
- 1. History: 3rd FIRF lower face
- 60/200/3.5/400 Neck centre, cheeks 2 passes (3 passes on lines) upper lip
- 60/200/2.0/400 chin
- **2. History:** 2nd FIRF lower face
- 50/200/3.5/400 Neck centre, cheeks 2 passes (3 passes on lines)
- 50/150/1.5/300 upper lip
- 50/200/2.0/400 chin
- Cold air
- tolerated well
- 3. History: 1st FIRF lower face
- intensitiy 50%
- Pulse duration 200ms
- Depth 3.0mm
- Delay 400ms
- just tolerated with cold air
- c/o numbing cream
- erythema

SECRET - Traumatic scars

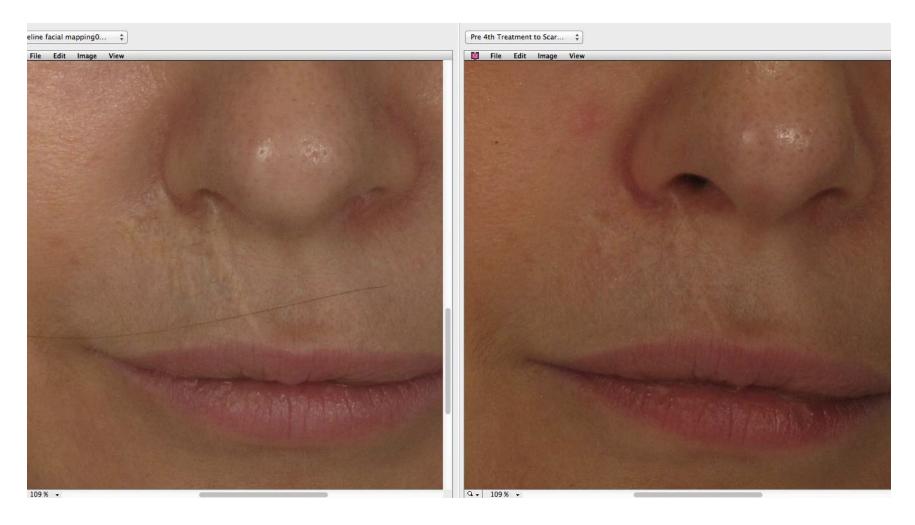

baseline 4mp3


- Presenting Problem: Traumatic scars
- 1. History: FIRF 75 shots
- 50%/3.5/200/400 20shots
- 50%/2.5/200/400 20 shots
- 50%/1.5/200/400 15 shots
- 2. **History:** secret previous was 70 shots
- today 3/2/1.5mm 75 shots
- 3. **History:** secret 40 shots (3 passes 50% overlap)
- 3.5/ 2.5/ 1.5 mm (300-600ms)
- 4. History: mca to r upper lip since 2 yrs old
- for sellas and secret x3 treatment

SECRET - Traumatic scars



baseline 4mp3


Acne scars Combo – 1550nm + SECRET + fillers

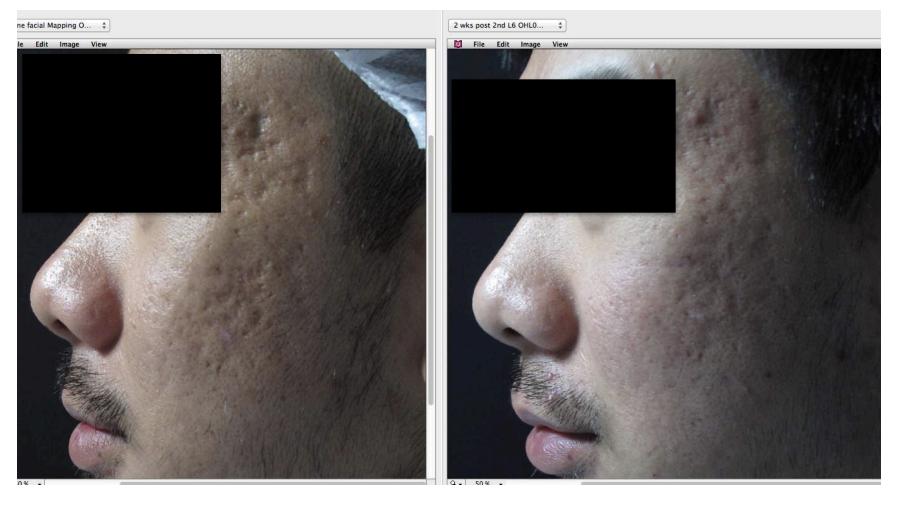

baseline 2wp2


- Combination Secret, 1550nm, fillers
- 1. History: fillers 0.25ml vol
- selas 2-3 passes 169ppa 90mj
- Secret:
- 60%/250/3.5/500 x2 then 150/3.5/300 x 1 pass cheeks 3 passes
- 60%/250/2.0/500 x2 then 150/2.0/300 x 1 pass temples 3 passes
- 60%/150/1.5/300 forehead 2 pass about 250 shots just tolerated with Siberian cooler
- about 230 shoots
- CA for 50 min

Combo – 1550nm + SECRET + fillers

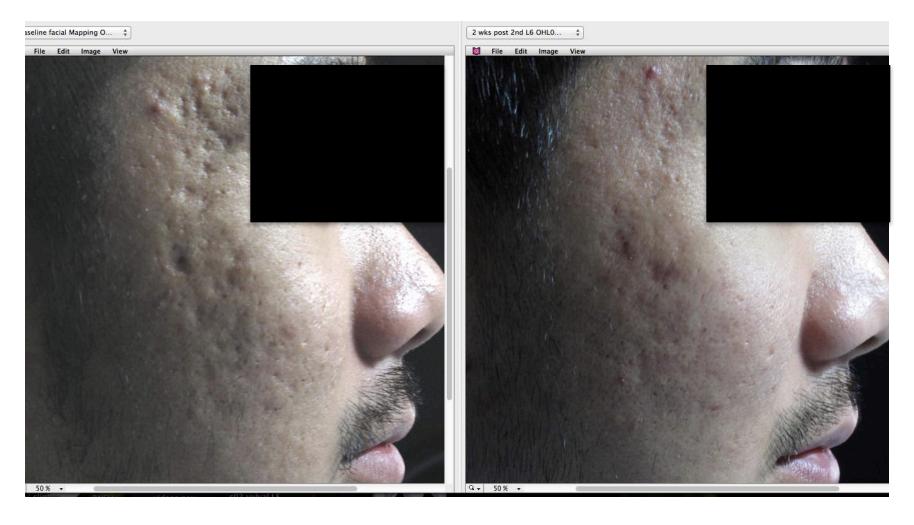

baseline 2wp2

Acne scar / melasma least likely to cause PIH

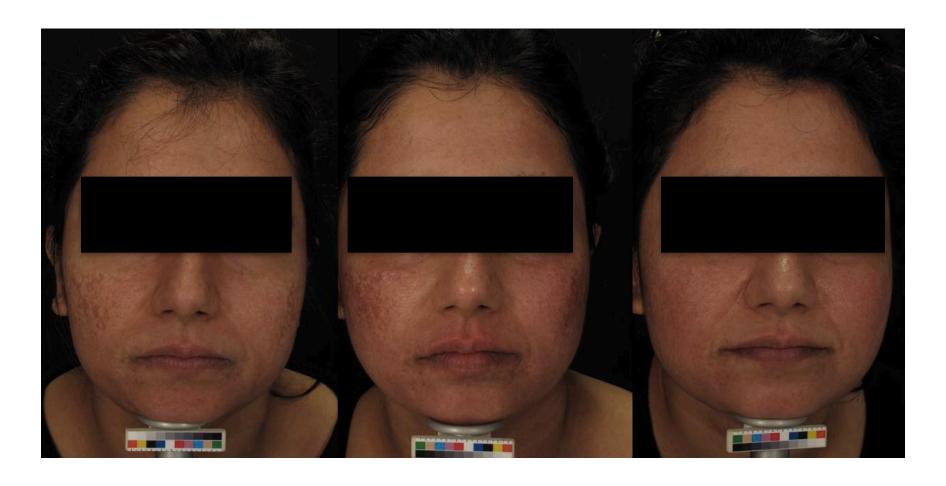

- Left Baseline: speckled dyschromia post QS for melasma
- Centre 3w post 2nd 1927: PIH
- Right 1mp1 Secret

- Secret (FIRF)
- 2.5mm / 250-500 50% 1 pass cheeks and NLF
- 1.5mm / 150-300 50% 2 pass cheeks and NLF